

Where comfort lives.

New to Respironics' Comfort Series™, the **ComfortClassic™** Nasal Mask raises the standards of comfort and convenience to deliver the right fit for your patients. From its new two-position gel spacer. To its reshaped membrane-thin silicone cushion. To an all-new "silent" exhalation port. It's more than a mask your patients live with. It's where comfort lives.

t h e c o m f o r t s e r i e s ™


ComfortClassicTM

The latest addition to The Comfort SeriesTM of masks delivers a more comfortable fit for your patients – and your business.

We designed the ComfortClassic to fit the same patient population with two mask sizes as previous masks did with three – so it's easier to stock. An integrated sizing gauge on the ComfortClassic's packaging helps in sizing and cuts down on mask waste by letting you determine the right mask size for your patient – before you open the bag.

For your patient, improved comfort starts at the point of contact. The ComfortClassic features a dual-thickness silicone cushion with an outer membrane that fills with air to create a comfortable seal. The thicker inner section provides stability.

A dual-position gel spacer molds to the wearer's forehead for even greater comfort and stability. ComfortClassic even features headgear attachment points positioned below the ear for more secure headgear placement. Combine these comfort features with an easy-to-maintain three-piece design, and you have a mask that may help to improve patient compliance. That's a feature anyone could be comfortable with.

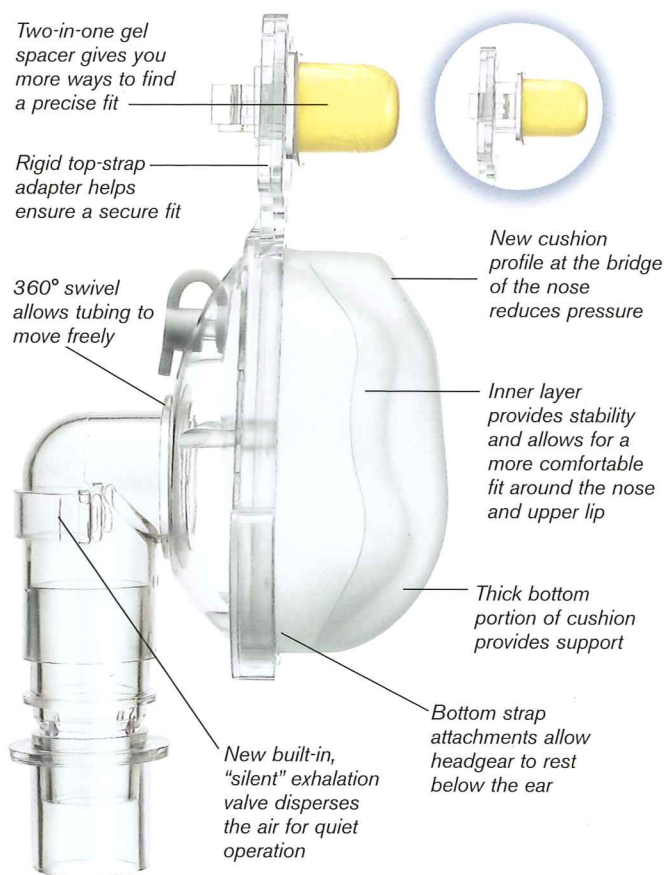

ordering info

	Part Number
SMALL WITH HEADGEAR	1007967
MEDIUM WITH HEADGEAR	1007968

replacement parts/accessories

DELUXE HEADGEAR ONLY	1002801
GEL FOREHEAD SPACER KIT	1002464
SPLIT WASHER KIT	1002720
EXHALATION ELBOW KIT	1004395
FOAM FOREHEAD SPACER, 5PK	1008651
SWIVEL	7041



ComfortClassic uses Respironics' one-size-fits-most **Deluxe Headgear**. This unique headgear with Velcro[®] tabs easily adjusts in eight places for a secure, custom fit. The Deluxe Headgear evenly distributes tension, adds stability and helps to prevent mask leaks. The ComfortClassic and Headgear disassemble easily for cleaning and replacement. The Headgear can be used with most Respironics' masks.




Customer Service: 1-800-345-6443 or 724-387-4000
 Respironics Europe: +33-(0)1-55-60-19-80
 Respironics Asia Pacific: +852-234-342-18


Latex Free • Patent Pending ComfortClassic and Comfort Series are trademarks of Respironics, Inc. Velcro[®] is a trademark of Velcro Industries ©2002 Respironics, Inc.
 CAUTION: U. S. federal law restricts this device to sale by or on the order of a physician.

1010280 SB 3/8/02


CARE
 Medical Limited
 Cardiac, Anaesthetic, Respiratory Equipment
 8H PIERMARK DRIVE, N.H.I.E, ALBANY
 TELEPHONE: +64 9 915-0800, FAX: +64 9 915-0801
 www.caremed.co.nz

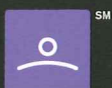

SLEEP
MANAGEMENT PROGRAM

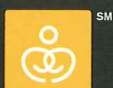

CHRONIC RESPIRATORY
MANAGEMENT PROGRAM

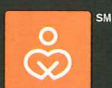

VENTILATION
MANAGEMENT PROGRAM



ASTHMA, ALLERGY & SINUSITIS
MANAGEMENT PROGRAM



INFANT
MANAGEMENT PROGRAM