


SleepSense®

Sensors that make sense

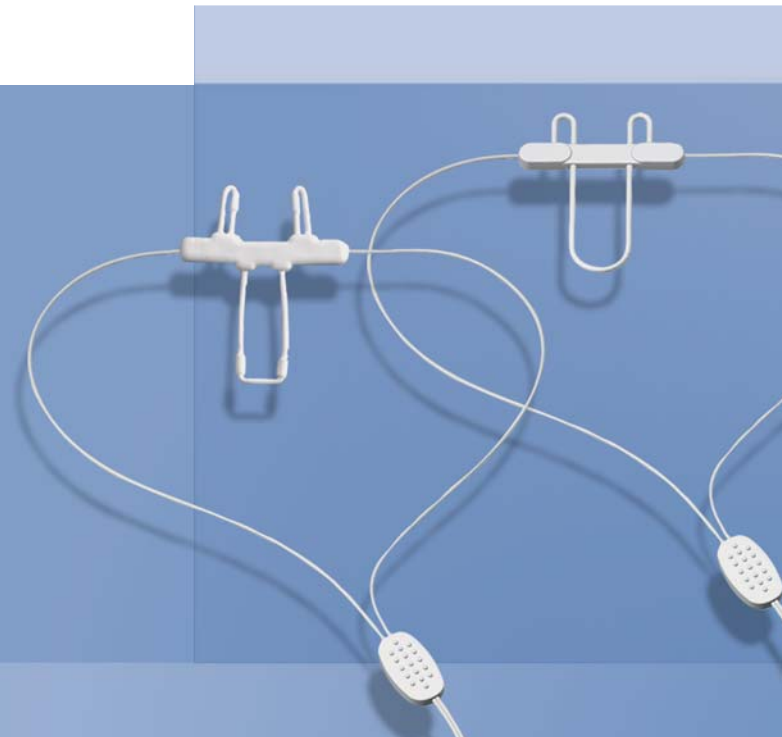

SleepSense® reusable flow sensors employ a high thermal mass design. This provides improved correlation with flow values, not just temperature.

These sensors are sturdy, stable, and will deliver a smooth flow signal throughout the night.

SleepSense® flow sensors share these advantages:


- ▶ Stiff prongs stay in place
- ▶ High thermal mass provides better correlation with flow, rather than temperature
- ▶ Small size and smooth contour for improved patient tolerance
- ▶ Optimized over 13 years of continuous development
- ▶ Robust and flexible, Kevlar® tinsel cables
- ▶ Sturdy mechanical design, 6 month warranty

SLP is certified to
ISO9001:2000 and
ISO 13485:2003 0483



**In the USA, FDA
restricts these devices
for sale by or on the
order of a physician**

SleepSense®

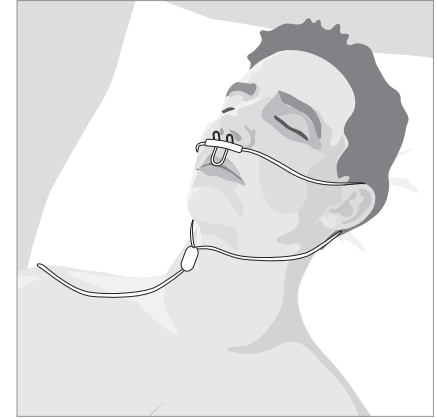

SleepSense® offers thermistor or thermocouple airflow sensors. Select the right technology and interface for your needs.

Thermistor sensors

These sensors provide the strongest and most accurate respiration airflow signal. They are especially robust, and easy to use and clean. They are powered by our interface box or your system.

1405	Thermistor flow sensor (P Series®)
1450-Kit	Thermistor flow sensor + interface box / safety pins
1451	Thermistor flow sensor (Suzanne®)
1453	Thermistor flow sensor (Edentrace®)
1459	Thermistor flow sensor (Alice 3®)
1460	Thermistor flow sensor (Alice 4®)
1461	Thermistor flow sensor (Alice 5®)



Note: 1450S (flow sensor only) is also available separately.

Thermocouple sensors

These sensors do not need an interface box, and are especially low noise due to very low internal resistance.

1401	Thermocouple flow sensor / safety pins
1402	Thermocouple flow sensor / 1/8" plug



Oral flow sensors

1413	Oral thermocouple flow sensor
1414	Oral thermistor flow sensor



Pediatric flow sensors

1400	Infant piezo thermal flow sensor / safety pins
------	--



For more information



www.sleepsense.com

S.L.P. Inc.

1121 E. Main St. #240, St. Charles, IL 60174, USA
 Tel: +1 630 513 7479, Fax: +1 630 513 8478
 E-mail: info@sleepsense.com
 Toll-free (U.S.): +1 888 757 7367