


# SleepSense®

Sensors that make sense

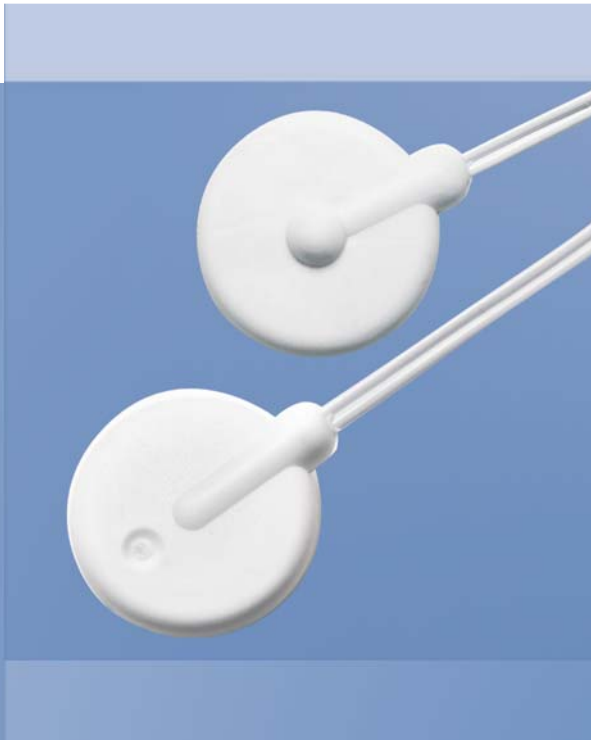

SleepSense® offers three sensing technologies for optimal matching with your clinical needs and system specifications.

SleepSense® snore sensors attach on the patient's neck to provide clear, artifact free qualitative indication of snoring and respiratory sounds.

## SleepSense® snore sensors share these advantages:

- ▶ Large sensitive area improves signal stability if sensor shifts during the night
- ▶ Soft, flat design has good seal against ambient noises
- ▶ Smooth contours and soft material protect the patient's skin
- ▶ High friction material decreases slide on oily skin
- ▶ Frequency response specifically tailored to snore monitoring
- ▶ Super strong and flexible, Kevlar® based cables
- ▶ Very durable, 6 month warranty


SLP is certified to  
ISO9001:2000 and  
ISO13485:2003



In the USA, FDA  
restricts these devices  
for sale by or on the  
order of a physician

## SleepSense®

SleepSense® snore sensors are available in three technologies. Select the right technology and interface for your needs.

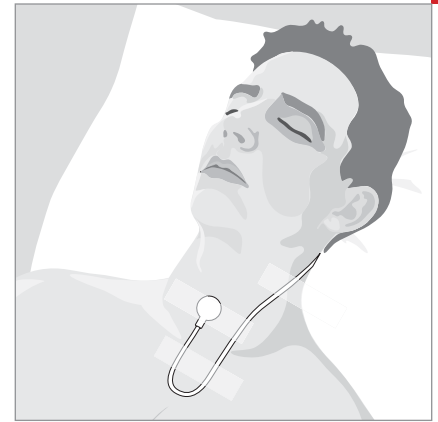


### Piezo sensors

Piezo sensors provide a signal proportional to the vibrations of the skin, and are immune to ambient noise. They are especially robust, and easy to use and clean.

1250	Piezo snore sensor / safety pins	
------	----------------------------------	-------------------------------------------------------------------------------------


### Dynamic sensors

Dynamic sensors provide a signal proportional to sound intensity, and have a very wide frequency response. They are especially suited for monitoring respiratory sounds in addition to snoring.

1200	Dynamic snore sensor / safety pins	
------	------------------------------------	---------------------------------------------------------------------------------------

### Electret sensors

Electret sensors are capacitive capsules powered by the system. They offer a wide frequency response and very high sensitivity. They are designed to match specific systems.

1205	Electret snore sensor (P Series®)	
1253	Electret snore sensor (EdenTrace®)	
1258	Electret snore sensor (Alice 3®)	
1261	Electret snore sensor (Alice 4®)	